


CHROMATOGRAPHIC SPECIALTIES INC.



Biphenyl Column, Acetonitrile Mobile Phase

Image 1

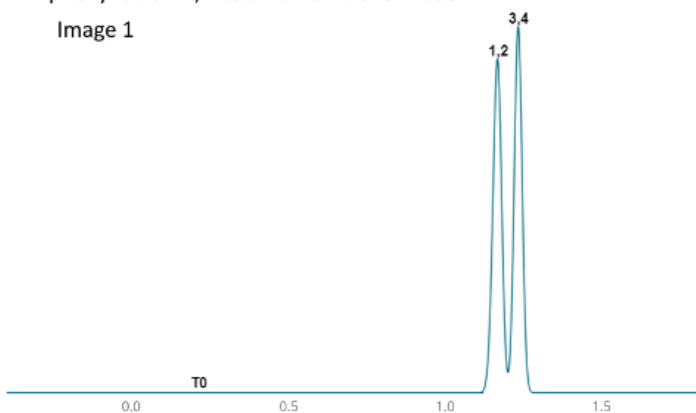

Biphenyl Column, Methanol Mobile Phase


Image 2

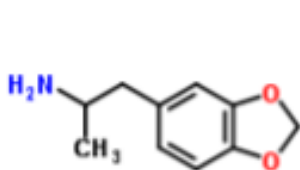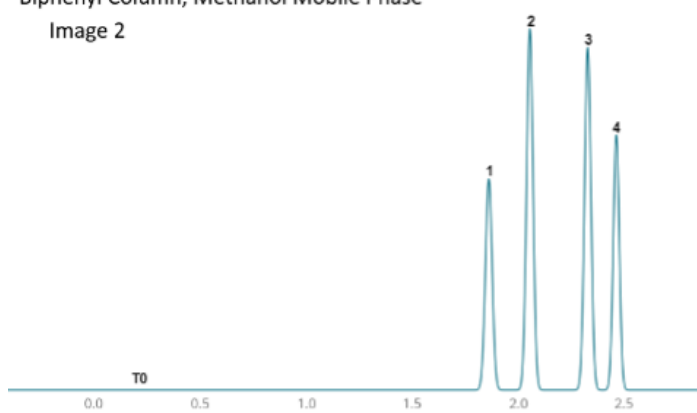


MDA (1)



MDMA (2)

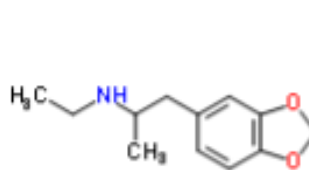

MDEA (3)

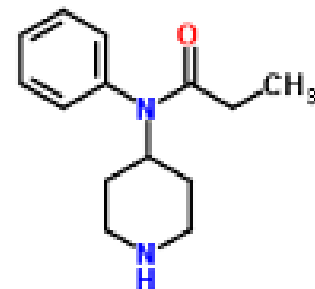

Norfentanyl (4)

Column Raptor Biphenyl (RK9309A5E)

Dimensions: 50 mm x 3 mm ID, 2.7 μ m

Mobile Phase:

A: Water, 0.1% Formic Acid

B: Acetonitrile (image 1) or methanol (image 2), 0.1% Formic Acid

Time (min)	Flow (mL/min)	%A	%B
0	0.8	96	4
5	0.8	0	100